

Mobile Sedation Carts

Configurable to fit your workflow

convergent[®] | innovative
medical

Efficient. Mobile. Customizable. Improve your bedside delivery workflow for patient sedation-analgesia with a dedicated cart that allows faster, more efficient care delivery, minimizes potential for error-causing work arounds and enhances patient safety and satisfaction. We developed this powered cart for a leading academic medical center to integrate med/narc storage, charting, patient monitoring and telemetry to perform bedside sedation from a compact, mobile platform used throughout their campus.




Applications:

- ✓ **Biopsy** liver, bone marrow
- ✓ **Imaging** X-ray radiography, ultrasound, CT, PET, nuclear medicine
- ✓ **Endoscopy** any procedure needed bedside for diagnosis or visualization
- ✓ **Pain management** sedation employed to relax a patient who is recovering from a painful event
- ✓ **Wound care** dressing changes
- ✓ For **pediatrics, impaired or disabled** patients who do not understand the situation

Mobile Sedation Carts

Configurable to fit your workflow

convergent[®] | innovative
medical

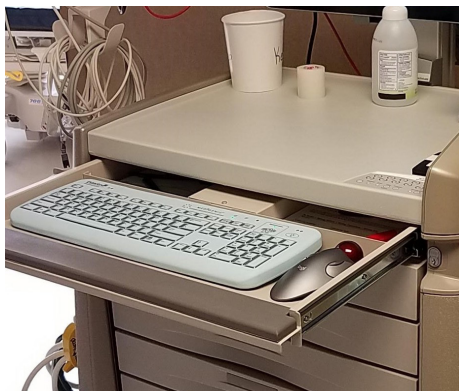

Features & Benefits

- ✓ **40AH Lithium Iron Phosphate Battery** – provides ample transit run time for any size campus, small or large.
- ✓ **Multi-LED Power Display** – Prominent display of cart battery charge status.
- ✓ **Keyboard drawer** – built-in with USB hub for quick, easy connection or replacement of keyboard, mouse.
- ✓ **Ample work surface** – includes a slide-out extension surface.
- ✓ **Enhanced Infection Control design** – minimizes surfaces & crevices to clean, places for fluids to pool.
- ✓ **Keyless, patented cart lock** – built-in SecureCatch™ feature, jam-proof design features and allows **selective drawer locking** capability for ‘always open’ drawers if desired
- ✓ **Badge-swipe access control** option available
- ✓ **Narcotics storage** options: PIN-access drawer or multiple ‘box’ options with BEST™ (replaceable) core lock
- ✓ **Steering Caster** – Easy-to-engage/disengage; it improves cart handling, minimizes steering effort when not at bedside—especially for traversing hallways in transit.
- ✓ **Central Display Support** – Accommodates up to 24” displays
- ✓ **CPU Compartment** – lockable, ventilated door; VESA and strap-mounting for variety of form factors.
- ✓ **Hospital-grade Power Cord & Multi-plug Power Strip**
- ✓ **VESA Mounts and Universal mounting channels** – Flexible mounting options for equipment, suction regulator/canister (note: hospital’s biomed responsible for fashioning a mounting plate for suction regulator and canister to Capsa Universal Mounting Channels)
- ✓ **Accessories** to support work flow: IV pole; sharps container; drawer dividers, bins, tubs, anesthesia trays, etc.
- ✓ **Enhanced Value and Total Cost of Ownership (TCO)** – features that keep the cart looking ‘new’ and functioning well beyond the life span of typical metal-box medical carts
- ✓ **Color-through polymer** body panels retain color if scratched; won’t dent
- ✓ **Molded-in handle** never loosens
- ✓ **Durable construction** – drawers operate smoothly for the life of the cart due to a welded steel tube bottom frame and jig-assembled aluminum upper frame to prevent ‘racking’
- ✓ **5-Year Warranty** (see full warranty for details)

